


Specifications

Product name	Mesh Nebulizer
Model	VP-M1
Power source	2 AA alkaline batteries(3V) Micro USB (5V)
Power consumption	Approximately 2 W
Nebulization rate	≥ 0.2 mL / min
Particle size	MMAD approximately 5 μm
Capacity of medicine cup	≤ 8 mL
External dimensions (W x H x D)	51 mm x 124 mm x 58 mm
Weight	Approximately 96 g (not including batteries)
Battery life	90 minutes with 2 AA alkaline batteries
Medicine cup life	180 hours
Working environment	Temperature: 10 to 40 °C Humidity: 30 to 85% RH
Accessories	Inhalation Mask (Children), Inhalation Mask (Adult), Mouthpiece, USB Cable
Pollution degree	Degree 2
Over voltage category	CATEGORY II
Altitude	≤ 2000 m

- Nebulization rate is lowered at low temperature.
- Specification and appearance may be changed without prior notice for improvement.



BMC's VP-M1 as a vibrating mesh nebulizer has an inner metal mesh vibrating in high-speed, makes the liquid medicine popup quickly when it passes through the thousands of laser-drilled holes. The small holes break medicine down into a very consistent particle size like PM2.5, forming the mist you inhale via the mask or mouthpiece. The mist can go into your lung directly to achieve a fast and effective treatment.

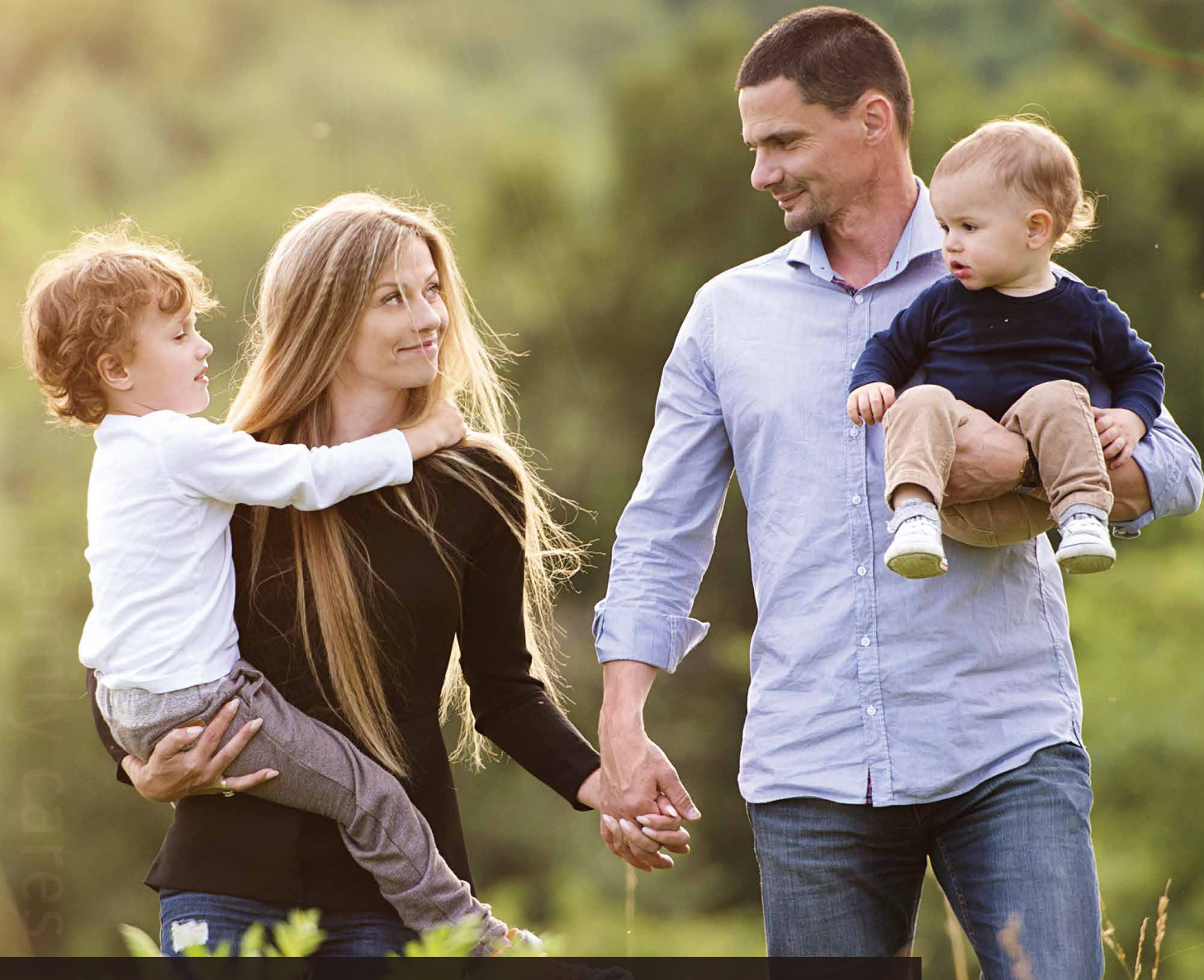

Mesh Nebulizer VP-M1



BMC Medical Co., Ltd.

Room 110, Tower A Fengyu Building, No.115 Fucheng Road, Haidian, 100036 Beijing, PEOPLE'S REPUBLIC OF CHINA

Tel: +86-10-5166 3880 / Fax: +86-10-5166 3880 ext. 810 / en.bmc-medical.com

ST-SC-20-02-11-23-001-Ver1.0

Mesh Nebulizer

Small Particles Size

VP-M1 works by a rapid vibrating mesh, the metal mesh products particles by very high-speed vibrating and makes the median particle size $\leq 5 \mu\text{m}$.

Medical Nebulizer

The maximum medication cap capacity can reach to 8 mL, meaning it can almost nebulize all the liquid medicine, more efficiency and professional.

Portable nebulization, held in one hand

Small and Light, not have to plug in power supply

The weight of the whole machine is about 96 g.
The dimension is 51 mm x 124 mm x 58 mm, VP-M1 is such small but smart. Only by two ordinary AA batteries for power supply, it can continue working for 90 minutes, which make the portable using be true.


Wisper-quite Operation

With a unique silent design, even if the baby is sleeping at night, VP-M1 can provide the nebulization therapy quietly.


Safe Material, no worries

The mask and mouthpiece part are made by food-grade materials. Suitable for all ages to use it without worries.

This portable Mesh Nebulizer offers a faster, more effective, more convenient medication administration in a quite comfortable way for both children and adult.

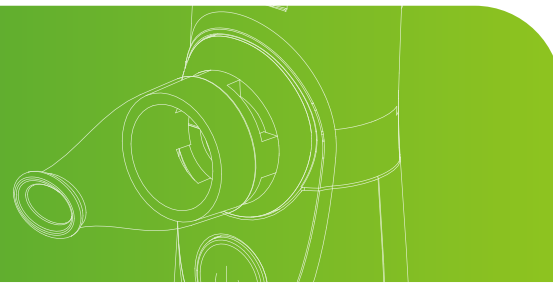

Name of Parts

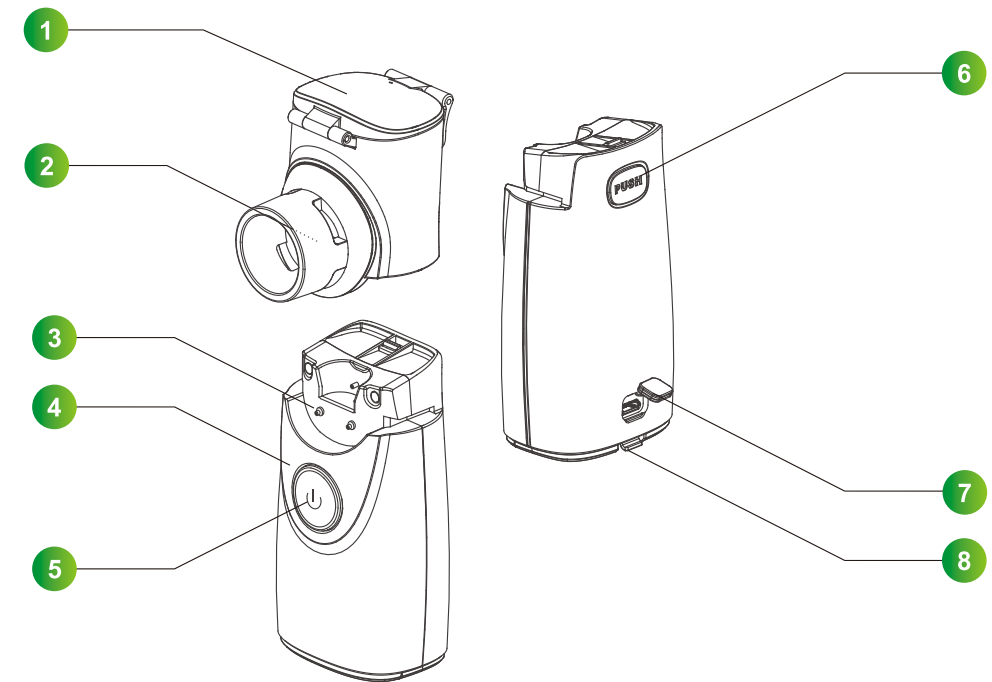

① Medicine Cup Cover

③ Electrode

⑤ Power Button

⑦ USB Cover


② Metal Mesh

④ Main Unit


⑥ Medicine Cup Button

⑧ Battery Cover

Contents

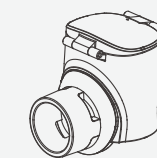

Medicine Cup

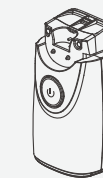

Main Unit




Mouthpiece




Inhalation Mask (Children)




Inhalation Mask (Adult)

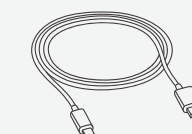

USB Cable

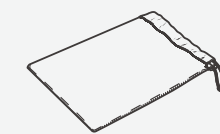

Storage Bag (Optional)

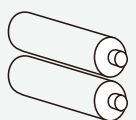

Battery